SEREFLO

synchroBreathe

Breath-Actuated Inhaler with dose counter

[S4] Reg. No. 42/21.5.4/0218

Sereflo 25/50 HFA

[S4] Reg. No. 42/21.5.4/0219

Sereflo 25/125 HFA

[S4] Reg. No. 42/21.5.4/0220

Sereflo 25/250 HFA

Instruction Leaflet

Please read this leaflet completely before using the inhaler.

Namibia: [NS2]

Sereflo 25/50 HFA :10/21.5.4/0507
Sereflo 25/125 HFA: 10/21.5.4/0505
Sereflo 25/250 HFA: 10/21.5.4/0506

Botswana: [S2]

Sereflo 25/50 HFA : BOT1101866
Sereflo 25/125 HFA: BOT1101867
Sereflo 25/250 HFA: BOT1101865

Cipla

CIPLA MEDPRO (PTY) LTD

Building 9, Parc du Cap,
Mispel Street, Bellville, 7530, RSA.


ABOUT THE

synchroBreathe

SynchroBreathe is an easy to use, breath-actuated inhaler. This means it releases a dose automatically when you breathe in. SynchroBreathe comes with a dose counter which helps you to check how many doses are left in the inhaler.

PARTS OF THE

synchroBreathe

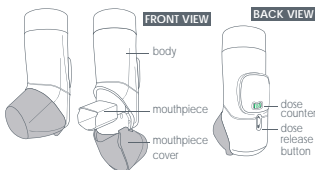HOW TO USE


synchroBreathe

I. Testing the SynchroBreathe:

If you are using the SynchroBreathe for the first time or it has not been used for a week or more, then shake it well and 'test fire' it by releasing two puffs into the air.

1. Shake and Open:

Hold the SynchroBreathe sideways, as shown in the picture below. Ensure that your fingers are away from the mouthpiece cover and not touching the dose release button at the back.

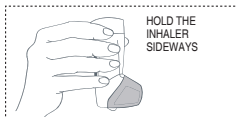
Shake your SynchroBreathe well, keeping the mouthpiece cover closed. Then Open your SynchroBreathe by placing your fingers on the grooves on the sides of the mouthpiece cover and gently pulling the cover down.


2. Fire a test dose:

Hold your SynchroBreathe upright and point the mouthpiece away from you. Release a puff by pushing down the dose release button which is located at the back of the SynchroBreathe just below the dose counter.

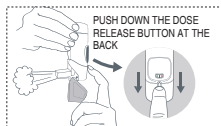3. Close:

Close the mouthpiece cover.

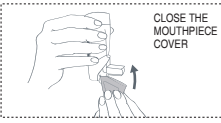To release another puff, repeat steps 1 to 3 above. Once 2 puffs have been released, your SynchroBreathe is now ready for use.

IMPORTANT:

- Test firing is ONLY needed if the inhaler is new or has not been used for more than one week.
- Always open the mouthpiece cover, gently.
- Always close the mouthpiece cover immediately after test firing. This prepares the SynchroBreathe for the next dose.
- Do not use the dose release button to take your medicine.

II. Using the SynchroBreathe


Step 1. Shake and Open:

- Hold the SynchroBreathe sideways, as shown in the picture below. Ensure that your fingers are away from the mouthpiece cover and not touching the dose release button at the back.
- Shake your SynchroBreathe well, keeping the mouthpiece cover closed. Then Open your SynchroBreathe by placing your fingers on the grooves on the sides of the mouthpiece cover and gently pulling the cover down.

Step 2. Breathe:


- While sitting or standing, hold the SynchroBreathe upright as shown.

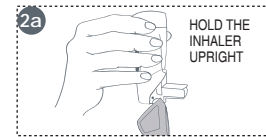
- Breathe out normally as far as is comfortable, through your mouth.

- Place the mouthpiece of the SynchroBreathe in your mouth; between your teeth and close your lips around it (do not bite). Start breathing in slowly and deeply through the mouthpiece. The SynchroBreathe will automatically release a dose when you begin to breathe in from the mouthpiece.

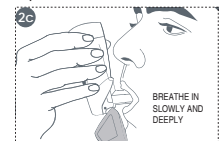
- Do not stop breathing in when you hear the click and whoosh, and feel the dose in your mouth. It is important to keep breathing in after the puff is released.
- Remove the SynchroBreathe from your mouth and hold your breath for 10 seconds, or for as long as is comfortable and then breathe out slowly.

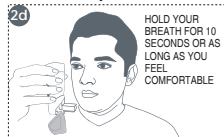Step 3. Close:

- After using your SynchroBreathe, hold it upright and close the mouthpiece cover immediately.

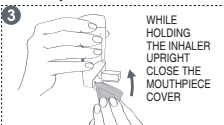Taking another puff:

Repeat steps 1-3 to take the next puff from your SynchroBreathe.

After you finish taking your puffs, rinse your mouth with water & spit it out. Do not swallow.

IMPORTANT:

- Always close the mouthpiece cover immediately after every puff. This prepares the SynchroBreathe for the next use.

HOW TO KNOW WHEN IT IS TIME TO GET A NEW

synchroBreathe

SynchroBreathe comes with a dose counter which is located at the back. The dose counter shows the number of puffs remaining. A new

SynchroBreathe will show the number 120 in the window, indicating the total number of doses present.

As you start using your SynchroBreathe, the dose counter will count down and indicate the number of puffs remaining.

Remember, the dose counter counts every dose but shows the number at intervals of 20. Thus, as the SynchroBreathe is used, the dose counter will count down and show the numbers 120, 100, 80, 60, 40, 20 and 0.

When there are 40 puffs remaining, the color of the number appearing in the window will change from green to red.

This indicates that only few doses are remaining in the SynchroBreathe.

Before the dose counter displays a '0' you should get a new Sereflo SynchroBreathe or ask your doctor if you need another one. When the dose counter displays '0', this means that there is no more medicine left.


Discard your SynchroBreathe once the number in the window displays a '0'. The SynchroBreathe may not feel empty and it may continue to operate, but you will not get the right amount of medicine, if it is used after it displays '0'.

CLEANING THE

synchroBreathe

- Clean the mouthpiece of your SynchroBreathe once a week.
- Open your SynchroBreathe by placing your fingers on the grooves on the sides of the mouthpiece cover and gently pulling the cover down.
- Wipe the mouthpiece with a clean dry cloth or tissue.
- Close the mouthpiece cover.

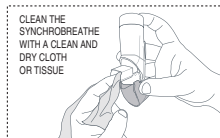IMPORTANT:

- Do not push a cloth or anything else into any part of the SynchroBreathe. This may cause damage to its operating parts.
- The SynchroBreathe should NEVER be washed. It should always be kept dry.

REMEMBER:

- Close the mouthpiece cover immediately after every inhalation.
- Do not open the mouthpiece cover with force nor rapidly as this may lead to the release of a dose.
- Don't open the mouthpiece cover by folding it down from the top as shown in the figure.

- Make sure that the mouthpiece cover is closed, when you shake the SynchroBreathe.
- If you are not sure how to use the SynchroBreathe, talk to your doctor.
- If the recommended dose does not provide the desired effect consult your doctor.
- You should never stop using your SynchroBreathe unless your doctor has told you to do so.
- Do not pull apart the SynchroBreathe.
- Take care not to drop the SynchroBreathe inhaler, as it may affect its functioning.
- Do not change/tamper with the number on the dose counter.
- Do not spray in the eyes.
- Do not puncture or burn the SynchroBreathe even when it is empty as it is pressurised.

Storing the SynchroBreathe:

Store in a cool dry place, at or below 25 °C. Avoid storage in direct sunlight and heat. The canister in SynchroBreathe is pressurised. The SynchroBreathe should not be broken, punctured or burnt, even when apparently empty.

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